

## EXHIBIT D

EIT UV Detector 365 nw  
Source Philips TL 29 D16/09n (Bulb)  
Ballast 120 vac with starter  
Detector 60 mm above Bulb  
Detector 30 mm above package  
Hole in pallet 20 mm  
Package 1 Day Acuvue


Readings in VDC	Empty package	Package with lens	Package with lens & packing solution
1	5.56	3.16	4.63
2	5.53	3.21	4.67
3	5.59	3.31	5.07
4	5.62	3.33	5.63
5	5.55	3.18	4.12
6	5.63	3.33	4.51
7	5.83	3.31	4.43
8	5.63	3.66	4.68
9	5.67	3.46	4.68
10	5.6	3.48	4.86
High	5.67	3.66	5.63
Low	5.53	3.16	4.12
Span	0.14	0.5	1.51
	Data A	Data C	Data E

EQ time  
Level of solution  
Package placement


Data A = Empty Package

Data C = Lens and Dry Package

Data E = Lens, Saline and Package

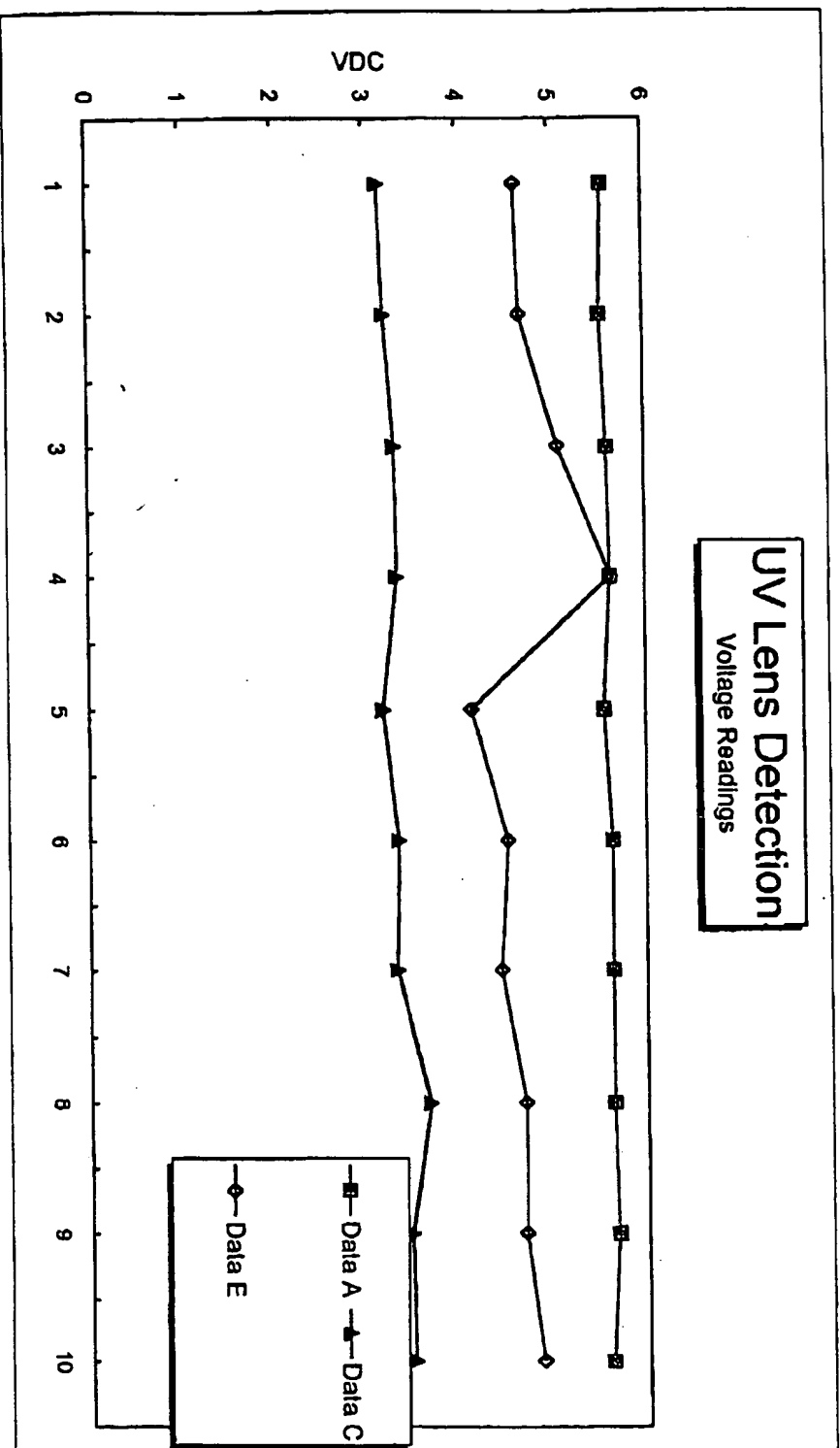

# UV-VIS

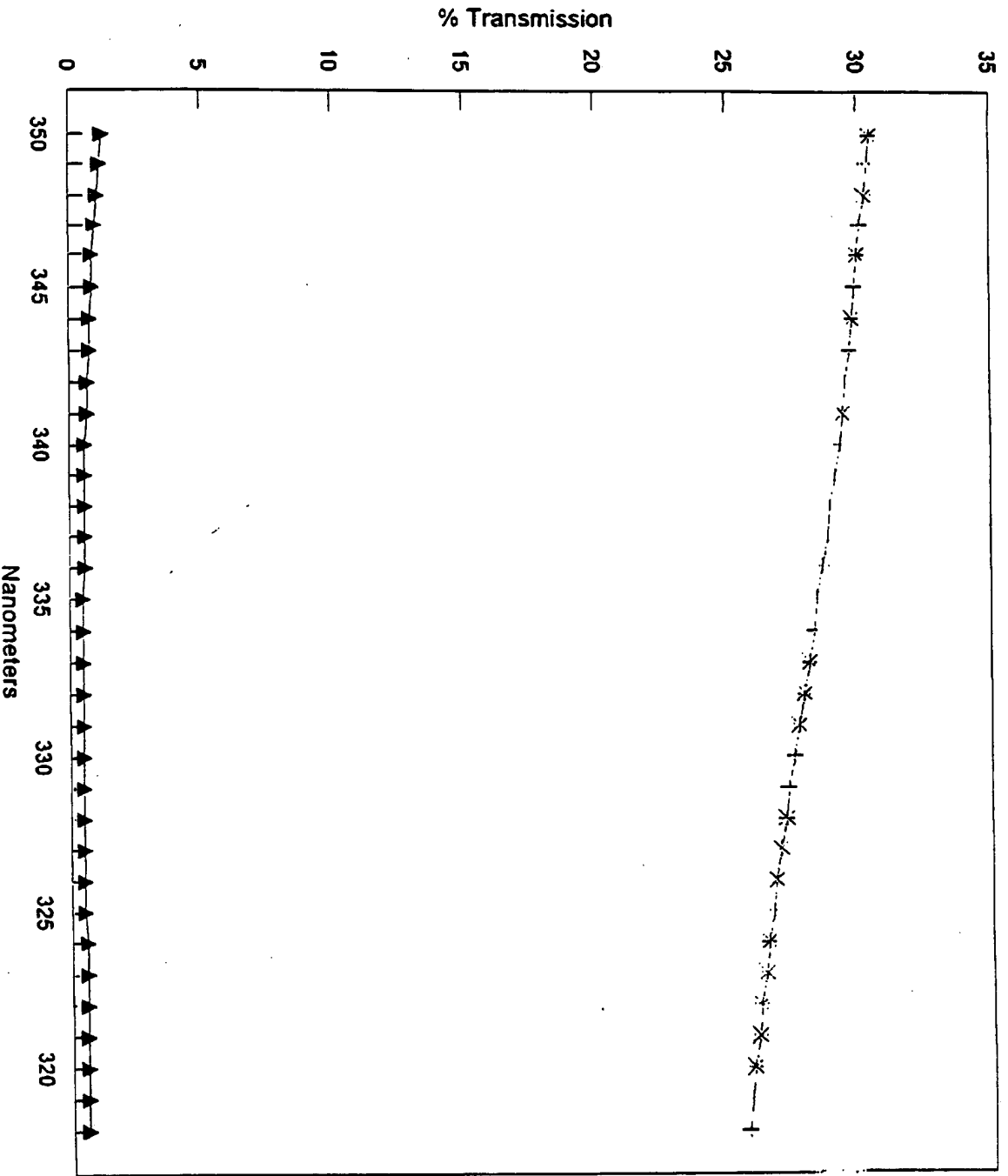

\* PACKAGE  
▲ LENS + PACKAGE



These pages are being submitted in  
Attorney Docket No.: VTN-0564  
Application No.: 09/819,074  
Office Action Dated: February 3, 2003

RECEIVED  
AUG - 6 2003  
TC 2800 MAIL ROOM

RECEIVED  
AUG - 7 2003  
TECHNOLOGY CENTER 2800